



Smithers AMS LLC
790 Main St
Wareham, MA 02571
License #: IL281355
Phone: 508-295-2550



Customer: The Healing Rose
Address:
23 Hale Street, Unit H
Newburyport MA 01950
Phone: 978-409-1091



Sample ID: 2051
Sample Name: Raw Cocoa 1x CBD Salve
Sample Type: Topical
Received Date: 08/31/2022
Report Date: 09/04/2022
Client Batch ID: Batch #1908
Client Sample Size: 63 g
Unit Weight (g): 63

Results Summary

Cannabinoids


09/04/2022

TESTED

Certificate of Analysis



Cannabinoid Distribution (%)

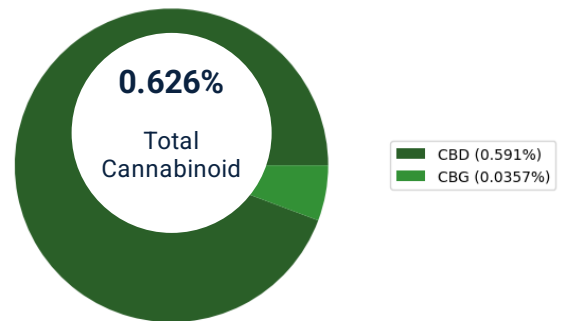

Disclaimer: The results presented within this report pertain only to the sample(s) indicated, as received by Smithers AMS LLC. Smithers AMS LLC is responsible for the information provided in this report, unless indicated that the customer provided the information. This report shall not be reproduced, except in full, without the written approval of Smithers AMS LLC. Action limits used to dictate pass/fail criteria are set by the Massachusetts Cannabis Control Commission.



Cannabinoids

Date Started: 09/01/2022
Instrument: UPLC-PDA

Date Completed: 09/04/2022
Method: TP-09.001.03 & TP-09.012.02

Analyst: Yüemin Wang Kovaleff
Deviation: N/A

Cannabinoid	Result (%)	Result (mg/unit)	LOQ (mg/unit)
CBD	0.591	372	0.00924
CBG	0.0357	22.5	0.00462
CBDVA	ND	ND	0.00369
CBDV	ND	ND	0.00369
CBDA	ND	ND	0.00924
CBGA	ND	ND	0.00462
Δ9-THCV	ND	ND	0.00369
Δ8-THCV	ND	ND	0.00462
THCVA	ND	ND	0.00369
CBN	ND	ND	0.00462
CBNA	ND	ND	0.00369
EXO-THC	ND	ND	0.00369
Δ9-THC	ND	ND	0.00924
Δ8-THC	ND	ND	0.00462
CBC	ND	ND	0.00369
CBL	ND	ND	0.00369
THCA-A	ND	ND	0.00924
CBCA	ND	ND	0.00369
Total THC	ND	ND	
Total CBD	0.591	372	
Total Cannabinoids	0.626	395	

Total THC = Δ9-THC + Δ8-THC + (THCA-A x 0.877)

Total CBD = CBD + (CBDA x 0.877)

Total Cannabinoids = Sum of all cannabinoids without THC/CBD corrections.

LOD = Limit of Detection; **LOQ** = Limit of Quantitation; **ND** = Not Detected.



Eric Steele, Lab Director
 09/04/2022

END OF REPORT